


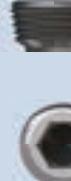
Product Catalogue

# Multi-Fix® SPACE SCREW



- Indication Prevention for Periosteum


This product is a disposable medical device and usage except for treatment purposes are prohibited.

SYSTEM	MINI PLATE	RADIUS PLATE	SMALL PLATE	LARGE PLATE TIBIA	LARGE PLATE FEMUR
모델명	SS2005CD	SS2505HD	SS3505HD	SS3500HD6	SS4500HD6
SCREW 호환직경	Ø 2.0	Ø 2.5	Ø 3.5	Ø 3.5	Ø 4.5
길이	0.5mm	0.5mm	0.5mm	0.0mm	0.0mm
체결부	CROSS	HEX 2.0	HEX 2.5	Hex.2.5	Hex.3.5
치료재료코드		BD1100AS	BD1100AS	BD1101AS	BD1101AS
체결방식					

0.0mm - 2.0mm 허가보유  
Material : Ti-6AL-4V-ELI ALLOY (ASTM F136)

### Prevention for Periosteum

There is a risk of blood supply disorder when the plate is closely stuck to the bone, which damages the periosteum. The purpose of Space screw is to prevent damage to the peritoneum by securing Plate and the bone space for the blood to flow smoothly.

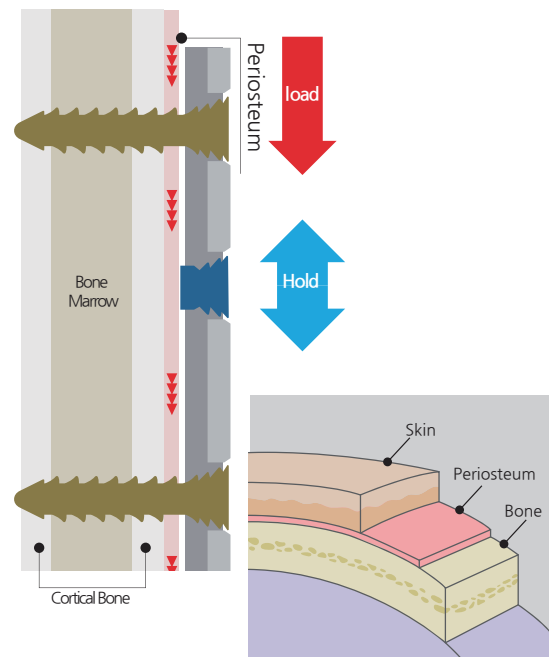

### Usage

1. Tighten the **PLATE durability by distributing the load that may occur in the empty HOLE area by fastening to the PLATE empty HOLE.**
2. By using Space screw to the empty HOLE of the Plate, the durability can be enhanced by distributing a load that may occur in the empty hole.
  - Supply and circulation of blood to fractured areas.
3. Prevent dislocation of bone graft.

### Advantages

1. Increase in stiffness, **strengthen load capability**
2. Can facilitate blood flow in the periosteum
3. Minimize unnecessary changes by being easily identifiable from other screws
4. Can be used as a customized product by having various specifications.
5. Case is easy to use and store.
6. Can use it smoothly and conveniently through pick-up driver.



Implant Case



Pick Up Driver